
SmartPump® Tourniquet System Rolling Stand

REF 5920-013-000

Instructions For Use

R_x ONLY

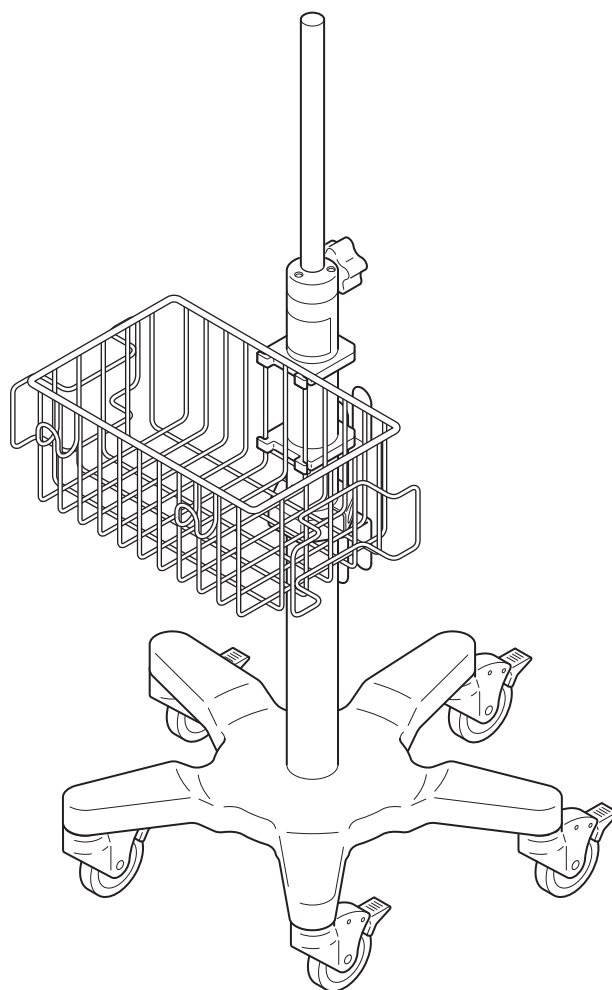ENGLISH (EN)
ESPAÑOL (ES)
DEUTSCH (DE)
FRANÇAIS (FR)
ITALIANO (IT)
NEDERLANDS (NL)
SVENSKA (SV)
DANSK (DA)
SUOMI (FI)
PORTUGUÊS (PT)
NORSK (NO)
POLSKI (PL)
ΕΛΛΗΝΙΚΑ (EL)
日本語 (JA)
中文 (ZH)
한국어 (KO)

Introduction

This *Instructions For Use* manual is the most comprehensive source of information for the use of your product. This manual may be used by in-service trainers and biomedical equipment technicians. Keep and consult this reference manual during the life of the product.

Conventions

The following conventions are used in this manual:

- A **WARNING** highlights a safety-related issue. ALWAYS comply with this information to prevent patient and/or healthcare staff injury.
- A **CAUTION** highlights a product reliability issue. ALWAYS comply with this information to prevent product damage.
- A **NOTE** supplements and/or clarifies procedural information.

Contact Information

For additional information, including safety information, or in-service training, contact your Stryker sales representative or call Stryker customer service at 1-800-253-3210. Outside the US, contact your nearest Stryker subsidiary.

Description

The Stryker Rolling Stand is designed to mount the Stryker SmartPump Tourniquet Pump in a safe manner. The stand allows for effective access to pump controls and relocation.

For Use With

This section identifies required components intended to be used with the equipment to obtain a safe combination. Component variations may exist.

WARNING: TIPPING HAZARD — DO NOT use the rolling stand REF 5920-013-000 with any other equipment, unless otherwise specified. ALWAYS use the rolling stand with the tourniquet pump. Failure to comply may result in patient and/or healthcare staff injury.

DESCRIPTION	REF
SmartPump Tourniquet System Single Channel Tourniquet Pump	5920-010-000
SmartPump Tourniquet System Dual Channel Tourniquet Pump	5920-011-000


User/Patient Safety

WARNINGS:

- Before using this equipment, or any component compatible with this equipment, read and understand the instructions for use. Pay particular attention to safety information. Become familiar with the equipment before use.
- Upon initial receipt, and before each use, inspect each component for damage. DO NOT use any equipment if damage is apparent.
- Use only Stryker-approved equipment, unless otherwise specified.
- DO NOT modify, service, or repair any equipment, unless otherwise specified. Call Stryker customer service. See the *Contact Information* section.
- Only individuals trained and experienced in the maintenance of this reusable medical device accessory should assemble and maintain this equipment.

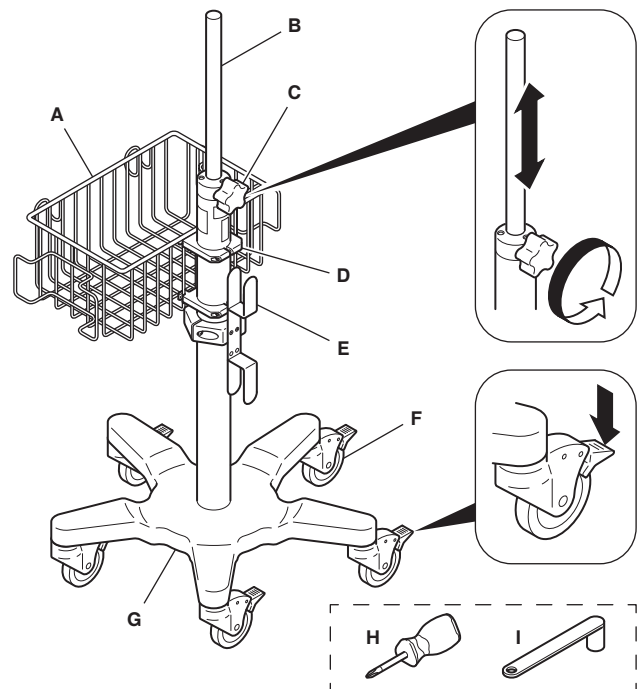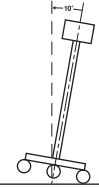
Features

A	Utility Basket
B	Post Assembly
C	Height Adjustment Knob
D	Basket Clips (2)
E	Power Cord Wrap with Clamp
F	Castors with Locks (5)
G	Base
H	Screwdriver
I	Wrench

Definitions

The symbols located on the equipment and/or labeling are defined in this section or in the *Symbol Definition Chart*. See the *Symbol Definition Chart* supplied with the equipment.

SYMBOL	DEFINITION
	To tighten the height adjustment knob and lock the post position
	To loosen the height adjustment knob and unlock the post position
	General warning sign
	 CAUTION: The loaded rolling stand must return to a vertical position when tipped 10 degrees. Read the instructions for use supplied with the rolling stand.

Assembly Instructions

NOTES:

- Refer to the illustrations and corresponding instructional steps.
- Make sure the following components are included before assembling the rolling stand: (1) utility basket, (2) basket clips, (1) power cord wrap with clamp, (1) post assembly, (1) base, (1) screwdriver, and (1) wrench. If any component is missing, call Stryker customer service. See the *Contact Information* section.

STEP	INSTRUCTION
1	Remove the hex head cap screw, lock washer, and flat washer from the post assembly. Use the wrench if necessary.
2	Fasten the base to the post assembly using the flat washer, lock washer, and hex head cap screw. Use the wrench to tighten the hex head cap screw. NOTE: Make sure the flat washer is against the base and the lock washer is between the flat washer and hex head cap screw.
3	Install and position each basket clip on the utility basket. Place the post and base assembly in an upright position. Remove the height adjustment knob. Install the utility basket and power cord wrap with clamp on the post assembly.
4	Install the height adjustment knob. Position the utility basket and power cord wrap on the post assembly as required. NOTE: Make sure the utility basket does not obstruct the labels on the post assembly. Use the screwdriver to tighten all the pan head machine screws to secure the basket clips and cord wrap clamp to the post assembly.

Inspection and Maintenance

See the *instructions for use* supplied with the SmartPump Tourniquet Pump.

Cleaning

See the *instructions for use* supplied with the SmartPump Tourniquet Pump.

Storage and Handling

See the *instructions for use* supplied with the SmartPump Tourniquet Pump.

Disposal/Recycle

See the *instructions for use* supplied with the SmartPump Tourniquet Pump.




Troubleshooting

See the *instructions for use* supplied with the SmartPump Tourniquet Pump.


Specifications

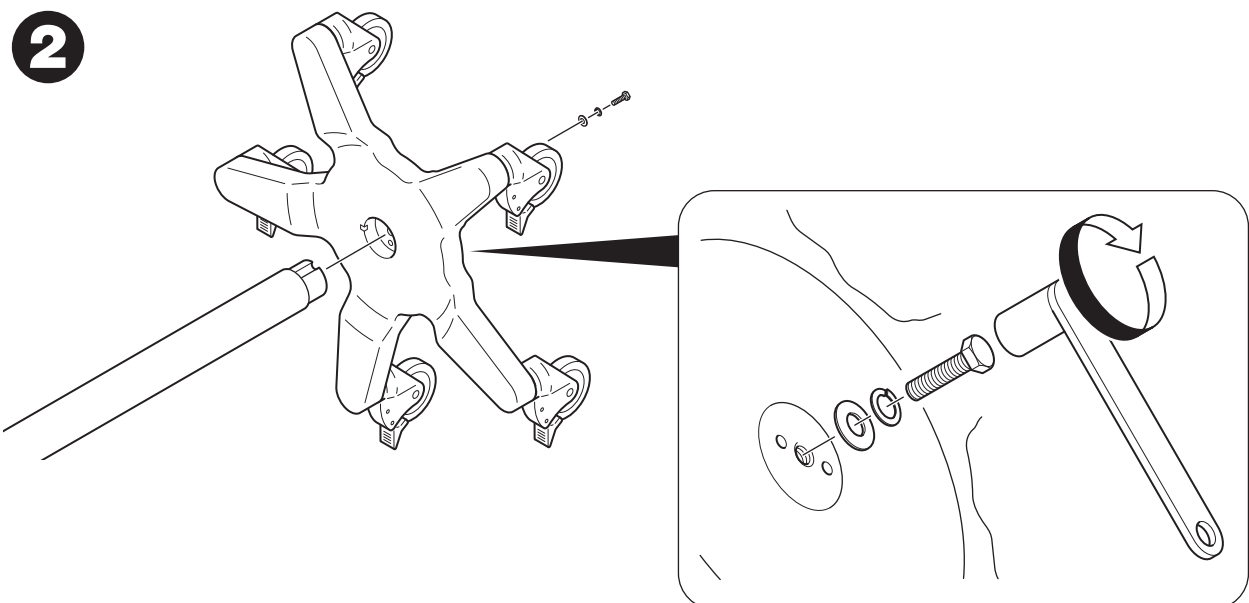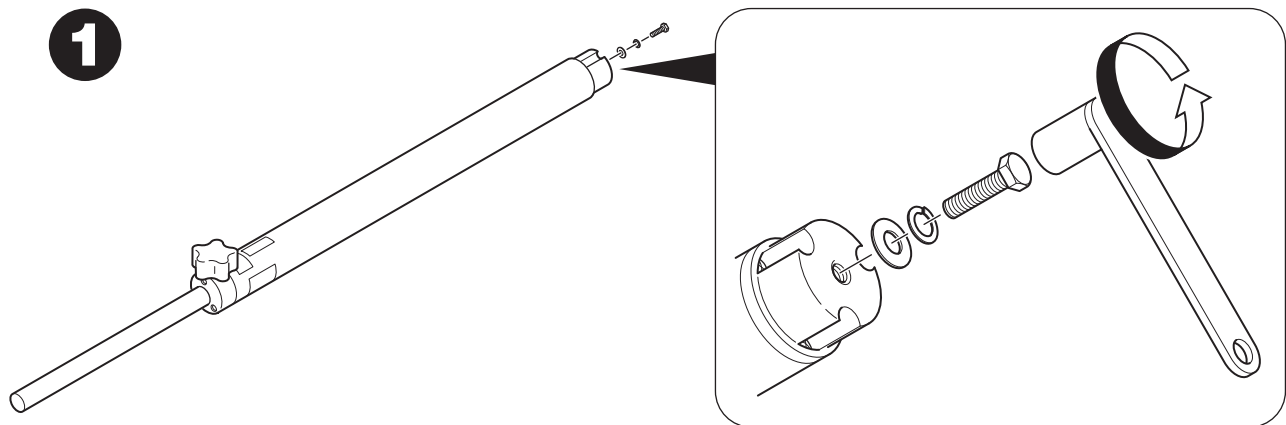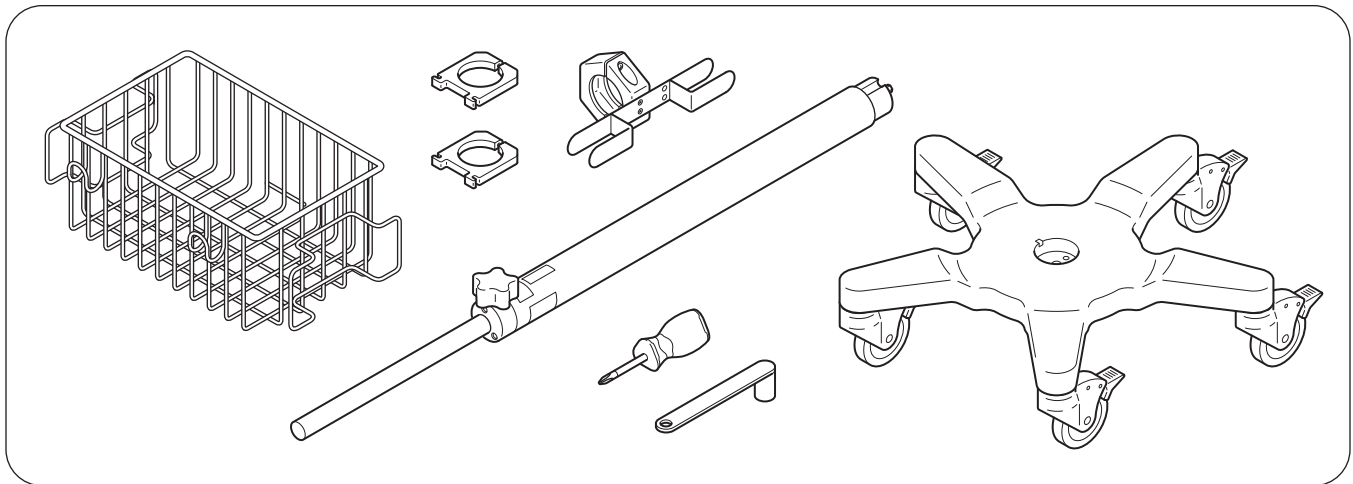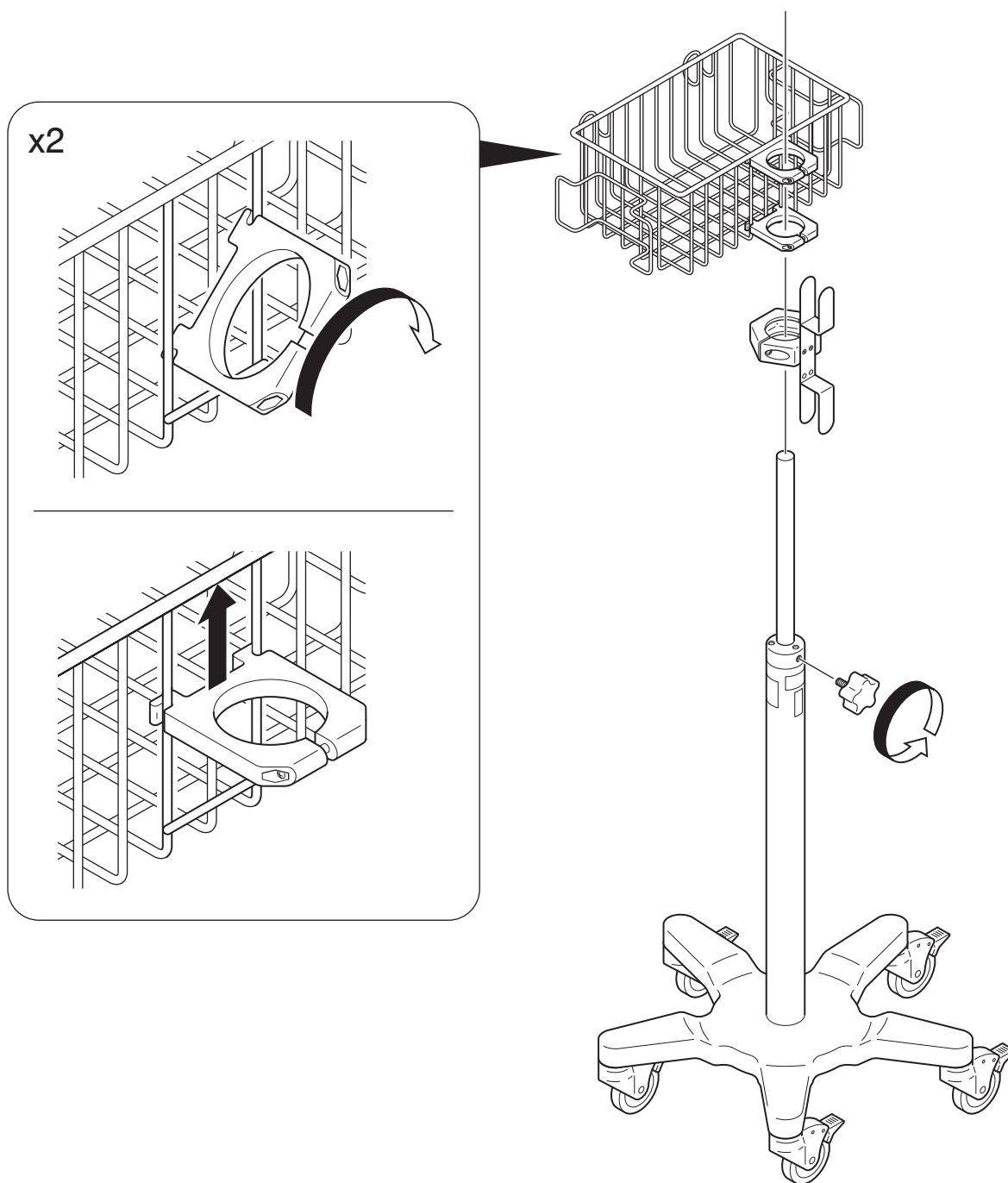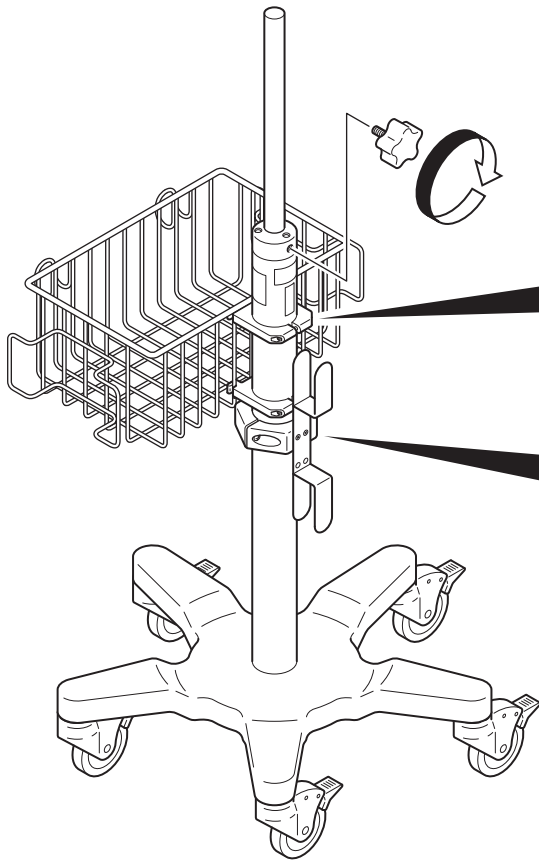
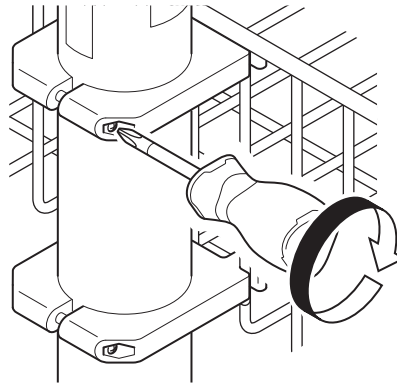
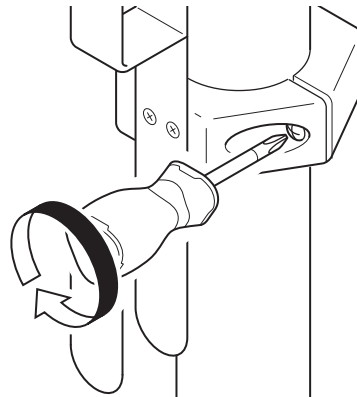
Model:	SmartPump Rolling Stand
REF:	5920-013-000
Mass:	11.2 kg [24.8 lb]; without any device; 4.5 kg [10 lb]; base counterweight
Dimensions	Height: 1.4 m [4.75 ft]; highest position; 1.1 m [3.8 ft]; lowest position
	Base Diameter: 0.5 m [1.75 ft]
Utility Basket Weight Limitation:	2.3 kg [5 lb]

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